

**Introducing... Starsys Reagent Shelving!**

**re-a-gent n.** A substance used in a chemical reaction to detect, analyze, or produce another substance.

Know as “reagent shelving” in the world of laboratory furniture, this new solid shelving option is now available for all sizes of Starsys Mobile WorkCenters! Starsys reagent shelving interfaces with slotted standards that are mounted within the vertical grooves of the same extruded aluminum channels that support overhead cabinets on Starsys Mobile WorkCenters. This new option is not field-retrofitable to existing units, it must be specified at time of order. Uprights are now available in three specific styles and *must* be specified for cabinets, slotted standards for reagent shelving or a combination of both.

**Construction:** The shelf assemblies are constructed of 14 gauge cold-rolled steel finished with a durable epoxy coating. Each shelf assembly consists of: one shelf, a pair of brackets and the associated hardware.

**Dimensional Information:**

Single - 21.5" (546mm)

Double - 41.3" (1049mm)

Triple - 61.1" (1552mm)

Single-Kneewell-Single – 73" (1854mm)

All shelves have a depth of 14" (355


Adjustable on 1" (25mm) centers

Minimum spacing between shelves is 5" (102mm)

**Load rating:** All four sizes of Starsys reagent shelving have been designed to meet the SEFA 8 shelf load rating standards of 50lbs per sq. ft. with a max load of 200lbs each. Total load for all Mobile WorkCenters is not recommended to exceed 900lbs, including unit weight.

**Assembly expectations:** The vertical supports and cross member are shipped individually packed and require assembly to the Starsys Mobile WorkCenter. (requires a T40 bit and .5" socket wrench) The Starsys Reagent shelf assemblies are individually packed. The side brackets require assembly to the main shelf body (requires a slotted screw driver). The shelf assembly can then be mounted at the desired location within slotted standards on one inch increments. When mounting 63" & 73" shelf assemblies it is recommended to have two people available for assembly.

**Lead Time** Consistent with the Starsys product system, reagent shelves are shipped within 4-6 weeks FOB Wilkes-Barre, PA USA.

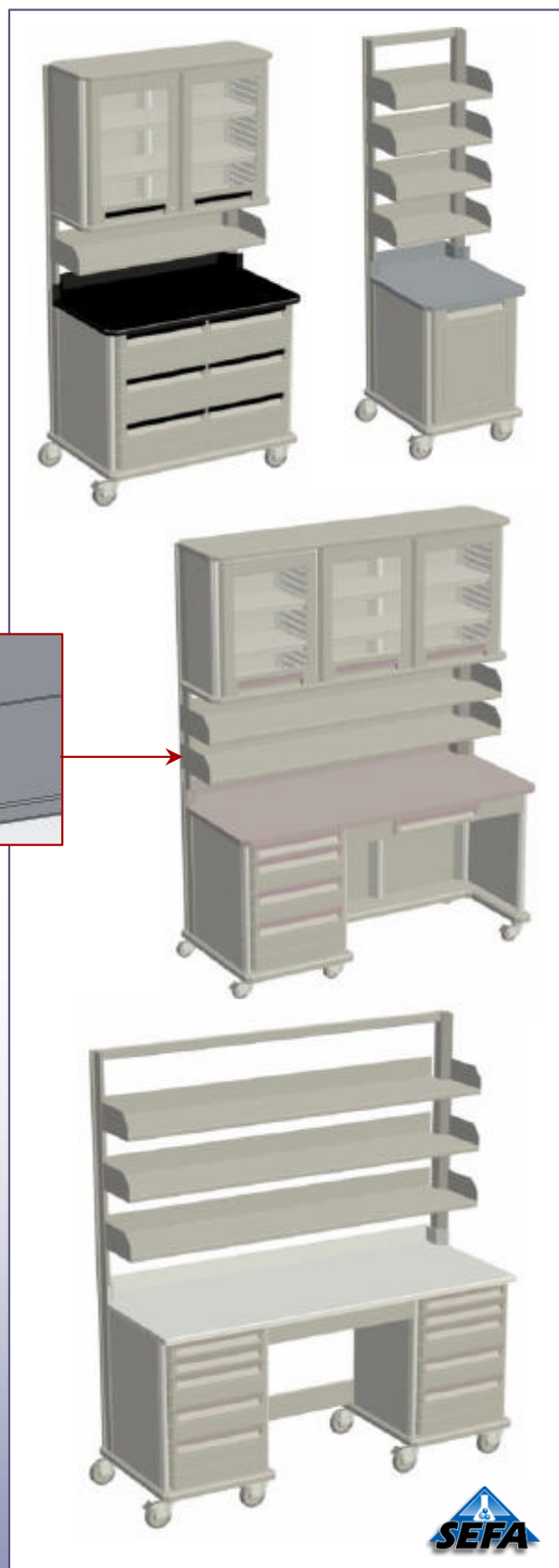



<b>Starsys Mobile WorkCenter Shelving &amp; Support Options</b>			
<b>Part Number</b>	<b>Description</b>	<b>Qty</b>	<b>Weight</b>
<b>Shelves</b>			
Reagent Shelf Assembles			
SXSMWSHF	Shelf Assembly - Single width	1	19
SXDMWSHF	Shelf Assembly - Double width	1	30
SXTMWSHF	Shelf Assembly - Triple width	1	42
SXKMWSHF	Shelf Assembly - Single-Kneewell-Single	1	47
<b>Lateral Supports</b>			
Lateral Support for MWC Shelving - 1 per MWC			
SXSMWSB	Cross member- Single width	1	4
SXDMWSB	Cross member- Double width	1	6
SXTMWSB	Cross member- Triple width	1	9
SXKMWSB	Cross member- Single-Kneewell-Single	1	10
<b>Lateral Supports</b>			
Lateral Supports for MWC Overheads - 1 pair per MWC			
SXSMWOHB	OH Structure Single Wide Mobile WC	1 pair	8
SXDMWOHB	OH Structure Double Wide Mobile WC	1 pair	12
SXTMWOHB	OH Structure Triple Wide Mobile WC	1 pair	17
SXKMWOHB	OH Structure Singlesw/Kneewell Mobile WC	1 pair	20
<b>MWC Uprights</b>			
MWC Uprights for Shelving ONLY			
SXMWUP24Q	24" base interior w/slotted inserts	1 pair	32
SXMWUP27Q	27" base interior w/slotted inserts	1 pair	32
SXMWUP30Q	30" base interior w/slotted inserts	1 pair	32
SXMWUP33Q	33" base interior w/slotted inserts	1 pair	32
<b>MWC Uprights</b>			
MWC Uprights for Overheads ONLY			
SXMWUP24-24	24" base interior with 24" OH interior	1 pair	28
SXMWUP24-27	24" base interior with 27" OH interior	1 pair	28
SXMWUP24-30	24" base interior with 30" OH interior	1 pair	28
SXMWUP24-33	24" base interior with 33" OH interior	1 pair	28
SXMWUP27-24	27" base interior with 24" OH interior	1 pair	28
SXMWUP27-27	27" base interior with 27" OH interior	1 pair	28
SXMWUP27-30	27" base interior with 30" OH interior	1 pair	28
SXMWUP30-24	30" base interior with 24" OH interior	1 pair	28
SXMWUP30-27	30" base interior with 27" OH interior	1 pair	28
SXMWUP33-24	33" base interior with 24" OH interior	1 pair	28
<b>MWC Uprights</b>			
MWC Uprights for Shelving AND Overheads			
SXMWUP24Q-24	24" base interior with 24" OH & slotted inserts	1 pair	32
SXMWUP24Q-27	24" base interior with 27" OH & slotted inserts	1 pair	32
SXMWUP24Q-30	24" base interior with 30" OH & slotted inserts	1 pair	32
SXMWUP24Q-33	24" base interior with 33" OH & slotted inserts	1 pair	32
SXMWUP27Q-24	27" base interior with 24" OH & slotted inserts	1 pair	32
SXMWUP27Q-27	27" base interior with 27" OH & slotted inserts	1 pair	32
SXMWUP27Q-30	27" base interior with 30" OH & slotted inserts	1 pair	32
SXMWUP30Q-24	30" base interior with 24" OH & slotted inserts	1 pair	32
SXMWUP30Q-27	30" base interior with 27" OH & slotted inserts	1 pair	32
SXMWUP33Q-24	33" base interior with 24" OH & slotted inserts	1 pair	32