

# ***Metro Lab Solutions***



# ***Presentation Overview***

---

- Storage Systems (Shelving)
- Lab Work Tables
- Lab Carts and Dollies (Polymer) : Utility, Glassware, Specialty
- Lab Carts Stainless Steel
- Autoclavable Animal Cage Racks and Casters
- Wall Shelving / Productivity Systems
- High-Density Storage Solutions
- Lab Furniture: Starsys™ Work Centers
- Gowning Room Solutions
- Mobile Computing Power System
- Metro Resources



# MetroMax i<sup>®</sup> Shelving System



- Corrosion Proof: Engineered polymer and Type 304 stainless steel construction. Lifetime warranty against rust and corrosion.
- Easy to Clean: Removable, lift-off shelf mats are a snap to clean, facilitating cGMP protocols.
- Antimicrobial Product Protection: Microban Antimicrobial product Protection - A Metro Exclusive – is built into all touch points.
- Durable: Strong, advanced polymer design supports up to 800 lbs. per shelf level.



# Super Adjustable™ Shelving System



- **Space Optimization:** Patented shelf release lever allows shelves to be quickly and easily adjusted in seconds.
- **Superior Air Circulation:** Open wire design minimizes dust accumulation and allows free circulation of air, greater visibility and light penetration.
- **Corrosion Proof:** Type 304 Stainless Steel finish, also available in chrome plated and corrosion resistant Metroseal 3™ epoxy.



# Super Erecta<sup>®</sup> Solid Shelving



- Choose between Flat and Louvered Embossed styles
- Type 304 stainless steel
- Available with stainless steel or epoxy coated cast corners
- Shelves feature a 1/8" raised "ships edge" on all four sides to contain spills



# ***Lab Work Tables: Stainless Steel***



- Stainless Steel construction
- Wide range of sizes, styles and accessories
- Mobile tables available, choose from numerous caster styles
- Table top selection includes Type 304 stainless steel and TRESPA Phenolic, available with or without backsplash. Perforated top tables available.



# *Lab Work Tables: MetroMax i*



- Advanced polymer construction provides superior corrosion resistance
- Wide range of sizes, styles and accessories
- Mobile tables available, choose from numerous caster styles
- Table top selection includes Type 304 stainless steel and TRESPA Phenolic, available with or without backsplash



# Lab Carts: MetroMax i™



Mobile Desk



Utility Cart



Chemical Transport

- Corrosion Resistant: Advanced polymers and Type 304 stainless steel construction
- Antimicrobial Product Protection: Microban Antimicrobial Product Protection is built-in
- Spill Retention: Solid bottom shelf contains spills and drips.



# ***Glassware Storage & Transport: MetroMax i***



- **Corrosion Resistant:** Advanced polymers and Type 304 stainless steel construction
- **Antimicrobial Product Protection:** Microban Antimicrobial Product Protection – A Metro Exclusive - is built into all pouch points.
- **Spill Retention:** Solid bottom shelf contains spills and drips.
- **Optional covers keep contents dust and dirt free**



# HPLC\* Carts: MetroMax i



- HPLC Carts free bench-top space increasing productivity
- Optional accessories include keyboard tray, AC power strip, Flat panel monitor mount, pull-out keyboard tray



\*HPLC – High-performance liquid chromatography

Architect and Designer Presentation



# Lab Utility Carts and Lab Dolly



- Corrosion Resistant: Polymer shelves are chemical and stain resistant, order units with stainless steel handles for a Corrosion Proof solution.
- Units available with stainless steel and chrome plated handles
- Microban Antimicrobial Product Protection: available with antimicrobial protection built-in
- Two and Three shelf models available.



# Lab Carts: Stainless Steel



Utility Cart

Wire and Solid Shelves  
Two and Three levels  
Type 304 Stainless or  
Chrome Plated



Security Cart

Secure storage for  
chemicals  
Type 304 Stainless



Autoclavable Cart

Two and Three levels  
Type 304 Stainless  
Wide choice of casters  
Autoclavable cart  
covers are available



# *Autoclavable Cage Racks*



- Stainless Steel construction
- Wide range of sizes, styles and accessories
- Choose from numerous caster styles

Autoclavable cart covers are available



# Autoclave Casters

- Hi-temp, autoclave casters (brake & swivel) for autoclave applications





Hi-temp Phenolic casters



Hi-temp Nylon casters



# SmartWall G3 – Productivity System

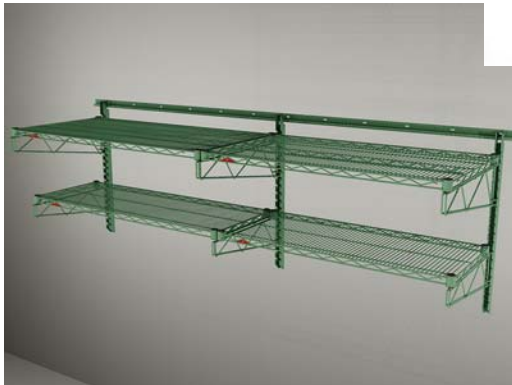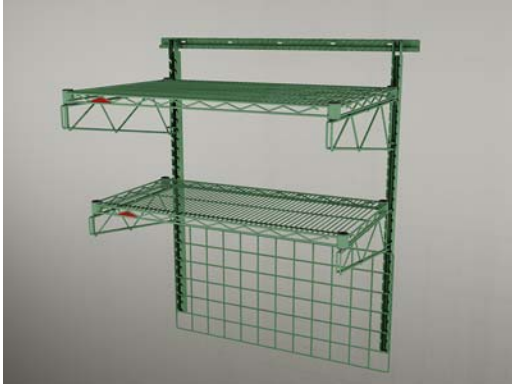

- Flexible wall storage
- Create stand-alone task stations with grids below shelves
- Bulk storage: robust weight capacities when installed into wood studs or equivalent structures
- Flexible runs of storage and task station configurations along the entire wall



# *Super Erecta Wall Shelving*



- Multiple tiers available
- 14", 18", 21", 24" depths
- Side by side units
- Stainless Steel, Chrome Plated and Metro Seal 3 available



**Load rating: 250 lb. per shelf not to exceed  
250 lbs. per unit.**

Architect and Designer Presentation



# Direct Mount Wall Shelving



- Typically used as single tier
- Side by side option available for SES
- Super Erecta style: Stainless Steel, Chrome Plated in 14", 18", 21", 24" deep
- MetroMax i style: 18" deep only

**Load rating: 200 lb. per shelf not to exceed 200 lbs. per unit.**



# Top Track™ High-Density System



- Increase storage capacity by up to 50%

- Mobile Units roll easily to create an “Active Aisle” between any two units

- Overhead track makes floor cleaning a snap



## Active Aisle



# qwikTRAK™ High-Density Storage



- Increase storage capacity by up to 50%

- Mobile Units roll easily to create an “Active Aisle” between any two units

- Heavy-duty weight capacity



**Active Aisle**



# ***Starsys Video***

## **Mobile, Modular, Metro Solutions That Move You!**



Double Click screen above to start Starsys Capabilities Video



# Starsys Laboratory Work Centers



- Available in Single, Double, and Triple widths
- Knee-well options available
- Configurable in multiple working heights: 33", 36", 39", 42"
- Multiple caster options available
- Countertop options include: Epoxy, Phenolic, Laminate, Solid Surface & Stainless Steel



# Starsys Laboratory Work Centers



- Available in a plethora of configurations
- Knee-well options available
- Configurable in multiple working heights: 33", 36", 39", 42"
- Multiple caster options available
- Countertop options include: Epoxy, Phenolic, Laminate, Solid Surface & Stainless



# *Starsys Laboratory Work Centers*



Architect and Designer Presentation

# Starsys Laboratory Work Centers



# *Starsys Laboratory Work Centers*




# *Starsys Laboratory Work Centers*



Architect and Designer Presentation

# Gowning Room Solutions

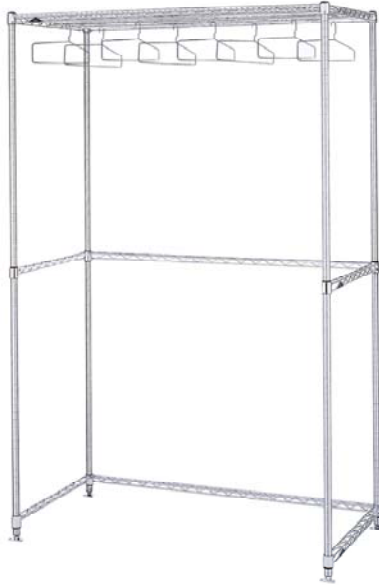


Garment Rack –  
Freestanding



Garment Rack – Wall  
Mounted




Gowning Bench –  
Floor Mountable or  
Freestanding

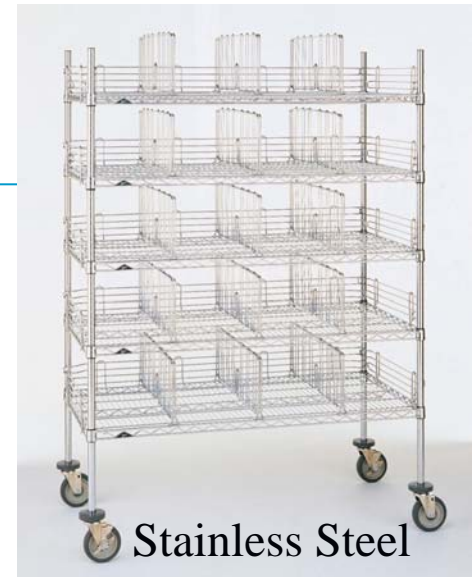

Stainless Steel

Gown / Bootie Storage

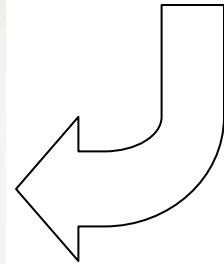

MetroMax I  
(Polymer)



# PowerSpot™ Mobile Computing Power System



- 6-8 hour run time: Integrated premium power system provides up to 150 watts of AC Power
- True Sine Wave Power:— Rapid rechargeable batteries convert to true sine wave AC output
- Corrosion Resistant: Mobile Power Units are constructed from advanced polymers and 304 stainless steel providing a lifetime of performance
- Maneuverable: Small footprint maximizes maneuverability
- Storage: Mobile Power Units have ample room for manuals, peripherals and other necessary items




# Metro Resources

- [www.metro.com](http://www.metro.com)
  - On-Line catalog
  - Space Planning Site
  - Web-based Product Configurator Tool (Password protected site)
    - Configure a unit in 5 minutes
- Architect and Designer Disc and Binder are packed with useful information
- Metro Representative
  - Metro representatives possess a vast array of application knowledge to customize a solution to meet your exact requirements
- Shelving Room Layout Tool
  - See next slide



NOTES:

- DIMENSIONS ROUNDED TO THE NEAREST 1/8"
- CUSTOMER IS RESPONSIBLE FOR SITE-Preparation TO ACCOMMODATE PRODUCT.

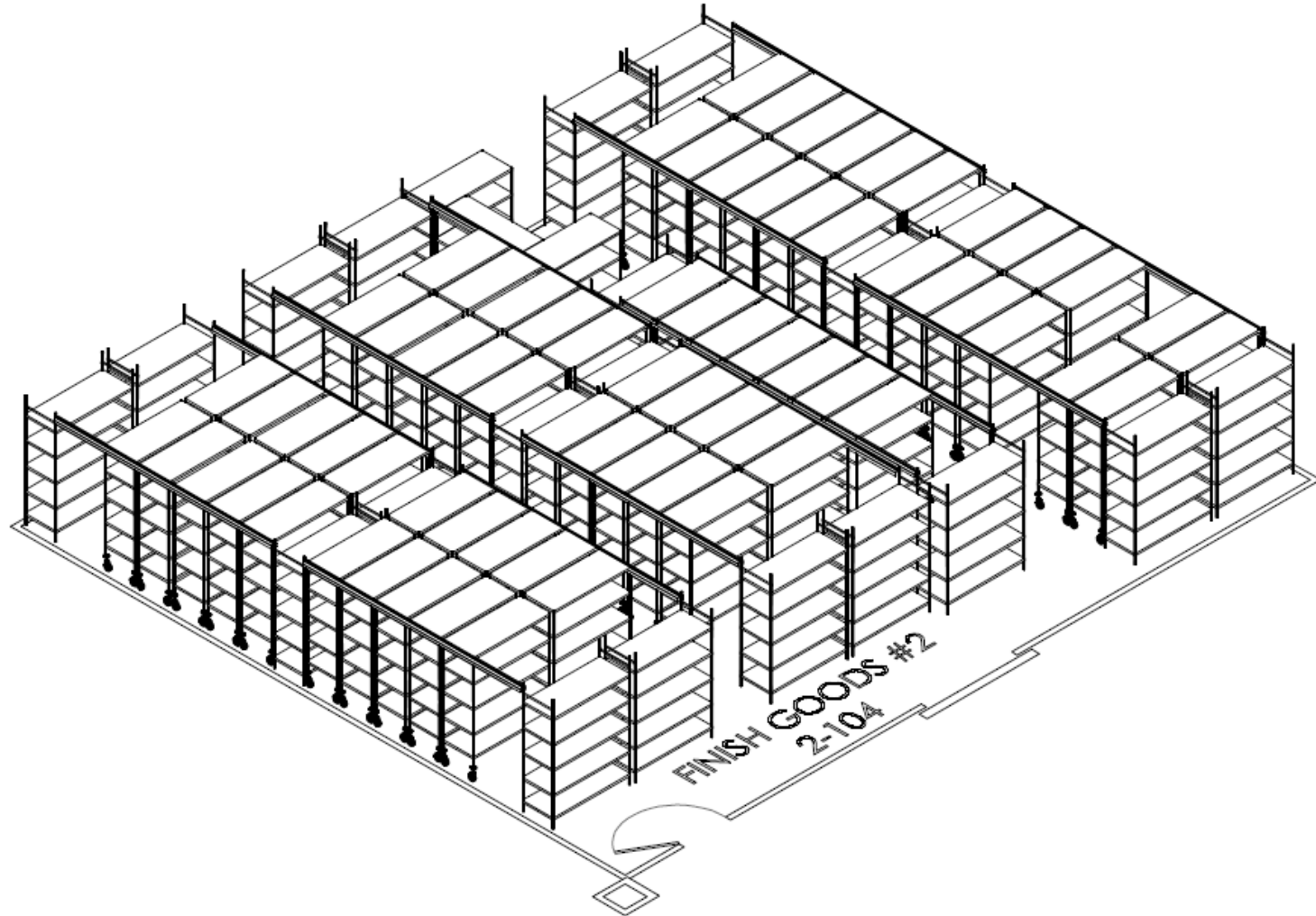

PLAN VIEW



InterMetro Industries Corp.  
651 N. Washington St.  
Wilkes-Barre, PA 18705  
PHONE: 570-706-4459  
FAX: 866-868-4120

CUSTOMER NAME: Canada Clean Room 12260  
PROJECT NAME: Abbott Healthcare Project  
DWG. TITLE: Finish Goods #2 2-104  
TERRITORY #: ORIGIN: J. Hallman

SCALE (ortho): 1/80  
SHEET#: 1 OF: 1  
DRAWN BY: FT  
REVISED:

FILE: 0108-Canada Clean Room  
QUOTE #:  
CREATION DATE: 1/10/08  
REVISION DATE: